

### **Attention! Important note!**

### **Read before replacing PCB KV-Control EZ130 (4512 178 00261)!**

This OPTIMUS RAD generator is a special generator. As opposed to the standard OPTIMUS RAD generator this special generator is not equipped with a RAD high-voltage tank, but it is operated with an R/F high-voltage tank and the associated expansion bellow.

This implies that PCB KV-Control EZ130 D800 (4512 178 00261) is not equipped with the KV-Control RAD-EPROM, but with a KV-Control R/F-EPROM (4512 113 26221).

In case PCB KV-Control EZ 130 (4512 178 00261) needs to be replaced **make sure to place the KV-Control R/F-EPROM (4512 113 26221) on EZ130 D800 of the new PCB KV-Control EZ130 (4512 178 00261).**

### Replacement Kit PCB KV-Control-3

4512 105 00051

### Level 0 Documentation

*This document and the information contained in it is proprietary and confidential information of Philips Medical Systems ("Philips") and may not be reproduced, copied in whole or in part, adapted, modified, disclosed to others, or disseminated without the prior written permission of the Philips Legal Department. This document is intended to be (a.) used by customers and is licensed to them as part of their Philips equipment purchase or (b.) used to meet regulatory commitments as required by the FDA under 21 CFR 1020.30 (and any amendments to it) and other local regulatory requirements. Use of this document by unauthorized persons is strictly prohibited.*

Printed in Hamburg, Germany

© 2005 Koninklijke Philips Electronics N.V.  
ALL RIGHTS RESERVED

**Proprietary Notice:**

*This document and the information contained in it is proprietary and confidential information of Philips Medical Systems ("Philips") and may not be reproduced, copied in whole or in part, adapted, modified, disclosed to others, or disseminated without the prior written permission of the Philips Legal Department. Use of this document and the information contained in it is strictly reserved for current Philips personnel and Philips customers who have a current and valid license from Philips for use by the customer's designated in-house service employee on equipment located at the customer's designated site. Use of this document by unauthorized persons is strictly prohibited. Report violation of these requirements to the Philips Legal Department. This document must be returned to Philips when the user is no longer licensed and in any event upon Philips' first written request.*

**Warranty Disclaimer Language:**

*Philips provides this DOCUMENT without warranty of any kind, implied or expressed, including, but not limited to, the implied warranties of merchantability and fitness for a particular purpose.*

**Limitations of Liability Language:**

*Philips has taken care to ensure the accuracy of this document. However, Philips assumes no liability for errors or omissions and reserves the right to make changes without further notice to any products herein to improve reliability, function, or design. Philips may make improvements or changes in the product(s) or program(s) described in this document at any time.*

## **SERVICE MANUAL – UNIT**

### **Replacement Kit PCB KV-Control-3**

**Author: Th. Frenschek**

File: Replacement kit PCB KV3\_988 00311AB

---

#### **List of pages and drawings (LOPAD)**

**Manual Order No: 4512 988 00311**  
**released: 06/2005**

1 ... 26      (05.1)

# Contents

<b>1.</b>	<b>INTRODUCTION AND TECHNICAL DATA .....</b>	<b>5</b>
1.1.	PURPOSE OF MANUAL .....	5
1.2.	ITEMS SUPPLIED .....	5
1.3.	COMPATIBILITY .....	5
1.4.	TOOLS/MATERIAL REQUIRED .....	5
1.5.	SAFETY INFORMATION .....	5
<b>2.</b>	<b>INSTALLATION .....</b>	<b>6</b>
2.1.	GENERAL .....	6
2.2.	PREPARATORY WORK .....	6
2.3.	CHECK OF THE GENERATOR TYPE .....	6
2.4.	INSTALLATION IN OPTIMUS RAD GENERATORS .....	7
2.4.1.	Installation of RAD software .....	7
2.4.2.	Adjustment of the factor for duty cycle .....	7
2.4.3.	Disposal .....	11
2.4.4.	Final test .....	11
2.5.	INSTALLATION IN OPTIMUS R/F OR OPTIMUS C GENERATORS .....	12
2.5.1.	Installation of the R/F software .....	12
2.5.2.	Adjustment of the factor for duty cycle .....	12
2.5.3.	Disposal .....	16
2.5.4.	Final test .....	16
2.6.	INSTALLATION IN OPTIMUS RAD FOR US ARMY APPLICATION GENERATORS .....	17
2.6.1.	Installation of the R/F software .....	17
2.6.2.	Adjustment of the factor for duty cycle .....	17
2.6.3.	Disposal .....	21
2.6.4.	Modification of the documentation .....	21
2.6.5.	Final test .....	21
<b>3.</b>	<b>TUBE CONDITIONING .....</b>	<b>22</b>
3.1.	GENERAL .....	22
3.2.	PRECONDITIONS / PROGRAM SETTINGS .....	22
3.3.	PROCEDURE .....	23

## 1. INTRODUCTION AND TECHNICAL DATA

### 1.1. PURPOSE OF MANUAL

When the service engineer orders an old PCB KV-Control-3 (4512 108 0908x), service logistics automatically sends a new PCB KV-Control-3 with the associated software PROM D800 4512 113 xxxxx according to the generator type OPTIMUS RAD, RAD for US Army application, R/F and C generators since the new PCB KV-Control-3 (4512 178 00261) is no longer compatible with the old KV-Control software.

### 1.2. ITEMS SUPPLIED

#### **4512 105 00051            Kit: "Replacement of PCB KV-Control-3"**

Consists of:

- |   |                |                         |               |
|---|----------------|-------------------------|---------------|
| - | 4512 178 00261 | Unit KV-Control-3       | (EZ 130)      |
| - | 4512 113 20141 | KV Control RAD software | (EZ 130 D800) |
| - | 4512 113 26221 | KV Control R/F software | (EZ 130 D800) |
| - | 4512 988 00311 | This manual             |               |

### 1.3. COMPATIBILITY

The kit "4512 105 00051" is compatible with PCB EZ130 "4512 108 0908x".

### 1.4. TOOLS/MATERIAL REQUIRED

- Standard tool set    TC 129

### 1.5. SAFETY INFORMATION

The general legal and factory safety recommendations for this X-ray equipment and the following recommendations must be strictly observed!

Start of installation, operation and maintenance work and especially electrical work must only be executed by trained and authorized persons. This equipment must only be serviced by properly educated service specialists who have received general and system-specific training as performed by Philips Medical Systems.



#### **WARNING**

*The system/component must be switched OFF during replacement work.  
Any X-ray unit produces ionizing radiation which may be harmful if not properly controlled. Therefore, it is recommended that this equipment be operated in accordance with the guidelines set down by the national council on radiation protection.*

---

## 2. INSTALLATION

### 2.1. GENERAL



#### WARNING

*When handling PCBs take all necessary electro static discharge (ESD) precautions!*

### 2.2. PREPARATORY WORK

- Turn OFF generator and ENF1.
- Remove PCB EZ 130 (4512 108 0908x).

### 2.3. CHECK OF THE GENERATOR TYPE

- Check on the mastercard, which type of OPTIMUS generator (RAD, RAD for US Army Application, R/F or C) is present on site:

OPTIMUS “**RAD generators**”  
contain  
**RAD HV tanks 9890 000 020xx.**

OPTIMUS “**R/F or C generators**”  
contain  
**R/F HV tanks 9890 000 027xx.**

“**OPTIMUS RAD for US Army  
Application generators**” are OPTIMUS  
RAD generators with KV-control  
**R/F software 4512 113 26212**  
(EZ130 D800) **and R/F HV tanks**  
**9890 000 00011** with expansion bellow  
(black).

9890-000-02183	BASE OPTIMUS	1
0001	MA-CONTROL	1
0002	CU BOOT	1
0003	CU INTERFACE CONTROL	1
0004	BASIS ZUBEHOER RAD - R/F	1
0005	Base Optimus	1
9890-000-02513	FIRMWARE OPTIMUS REL.3.6	1
0001	CENTRAL UNIT REL.3.6	1
9890-000-02051	HE ERZEUGER 1T 65/80/100 KW	1
0001	HE-H 1 ARB. XRG 90 RAD	1
9890-000-02811	AEC AUTOMATIC EXPOSURE CONTR.	1
0001	SW-OPTION ACCORD. TO CONF. KEY	1
9890-000-02821	APR/F ANATOMIC PROGR.R/F	1
0001	APR ACCORD. TO CUSTOMER ORDER	1
9890-000-02091	65/80 KW RD EXTENSION 480 V	1
0001	KV-CONTROL RD	1
0002	AUSBAUSATZ 65/80 KW RAD, Tr.	1
9890-000-02561	FLAENCHENDOSIS KALKULATOR OM	1
0001	SW-OPTION ACCORD. TO CONF. KEY	1
9890-000-02211	VARIOFOCUS OPTIMUS	1
0001	SW-OPTION ACCORD. TO CONF. KEY	1
9890-000-02332	1 ADAPTER AMPL.KABEL OPT.	1
0001	AMPLIMATKABEL-ADAPTER	1
9890-000-02472	SURGE ADAPTER OPT	1

High-voltage tank 9890 000 0xxxx

- With OPTIMUS RAD generators refer to 2.4. (page 7).
- With OPTIMUS R/F or OPTIMUS C generators refer to 2.5. (page 12).
- With OPTIMUS RAD for US Army Application generators refer to 2.6. (page 17).

## 2.4. INSTALLATION IN OPTIMUS RAD GENERATORS

### 2.4.1. Installation of RAD software

- Install the KV-Control RAD software 4512 113 20141 on the new PCB EZ 130 (4512 178 00261) in socket D800.
- Install the new PCB EZ130 (4512 178 00261) in the generator.

### 2.4.2. Adjustment of the factor for duty cycle

#### 2.4.2.1. General information


The actual value of the set kV must be attained at least after 2ms. During the kV rise phase there must be neither kV break-in nor a kV overshoot.

The factor duty cycle is based on an adapted tube and determines at local mains voltage and mains resistance conditions:

- kV rise phase
- kV behavior during the exposure in falling load technique.

The factor duty cycle is stored in the memory of PCB CU EZ139.

Refer to figure 1:

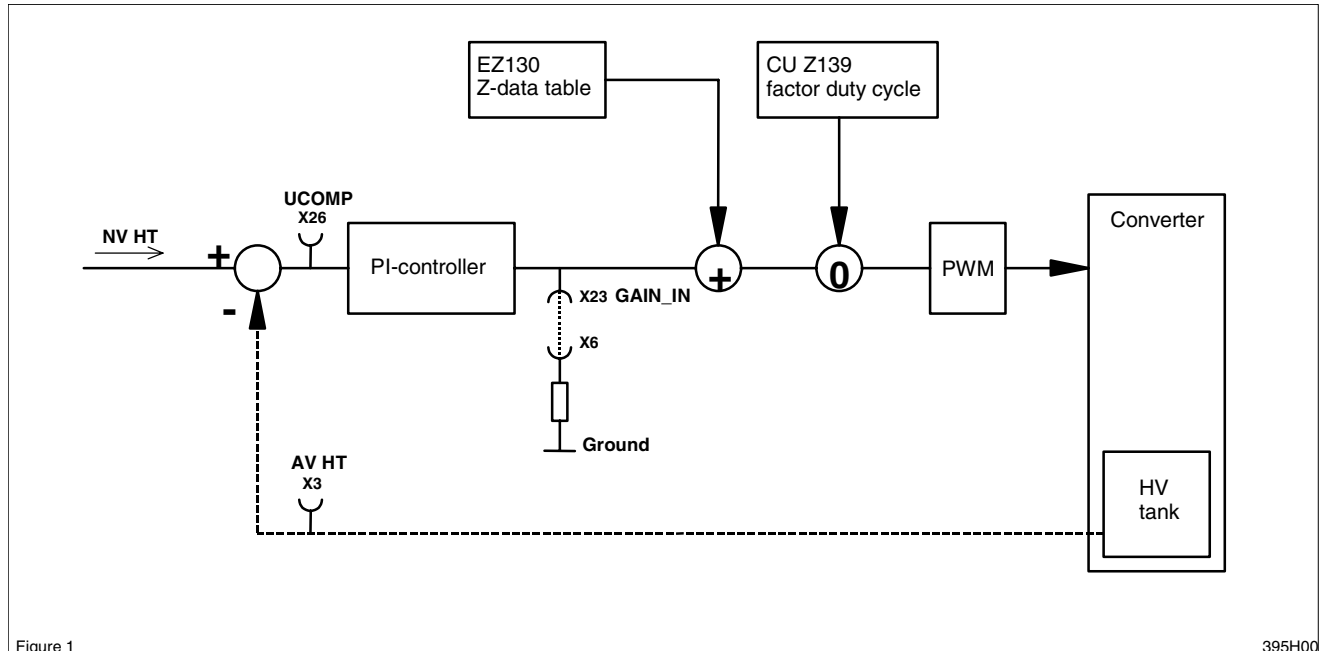


During alignment this factor duty cycle must be entered via AGenT. The influence of this factor as a correction value for the Z-data table is monitored as the  $U_{COMP}$  signal, since the PI-controller is deactivated by the grounded  $GN\_IN$  signal.



### 2.4.2.2. Connecting and setting the scope

For connections see figure 2:

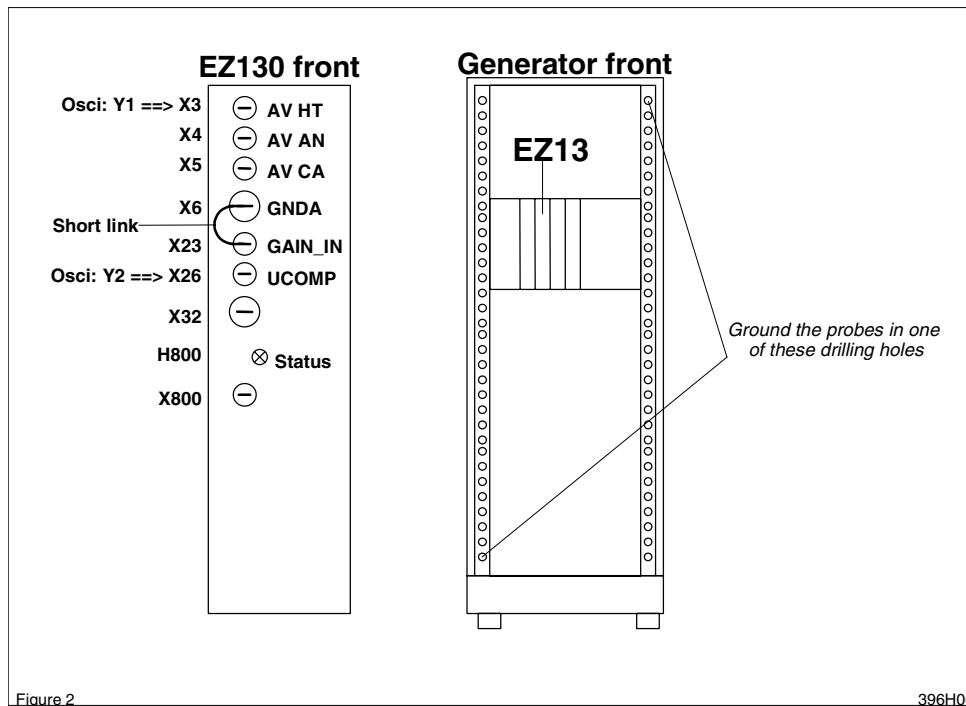

Channel 1 = EZ130 X3 ---> AV HT ---> 20kV/V ---> 1V/div --> Zero-line at bottom of screen  
 Probe GND = one of the drilling holes at the front cabinet chassis

Channel 2 = EZ130 X26 --->  $U_{COMP}$  ---> 1V/div ---> Zero-line 2 div from bottom of screen  
 Probe GND = one of the drilling holes at the front cabinet chassis

Trigger = external (preferred) ---> CTRL\_X\_C/ ---> backpanel EZX74 / negative slope  
 or = internal channel 1 ---> AV HT ---> EZ130 X3 / positive slope at +3V  
 Probe GND = one of the drilling holes at the front cabinet chassis

Time base = 5 or 10ms/div ---> trigger delay -1div



#### NOTE

*A digital scope should not have any other ground connection than the ground of the three probes at the drilling holes at the front generator chassis.*

*A mains-driven scope must be isolated from earth potential, otherwise it might display artefacts.*

#### 2.4.2.3. Deactivating the kV controller

- Connect EZ130 X23 GAIN\_IN and X6 on board with a short link (use a short wire).



#### CAUTION

*This alignment requires exposures with high kV.  
Be sure the tube has been warmed up before.*

---

#### 2.4.2.4. Setting of exposure data

- Switch ENF1 and the generator ON.
- Set 125kV.



#### NOTE

*Any tube arcing during this adjustment requires the execution of the tube conditioning next as described in chapter 3 (see page 22).*

*Disconnect the short link between X23 and X6.*


*Start over this adjustment from chapter 2.4.2.3. onwards if the tube conditioning was successful.*

---

- Set kV and mA values according to the programmed tube limits:  
**125kV:** 100mA at kV\_3 (50kW).
- Set the exposure time: 40ms.

#### 2.4.2.5. Adjustment of the factor for duty cycle

- Adjust the factor duty cycle via service software AGenT by measuring  $U_{COMP}$  with the scope.
- Connect the service PC and start AGenT:  
Select menu:  
*Program / FU KV IGBT Pulse Width correction.*
- Set the starting value factor duty cycle to **1.00**:

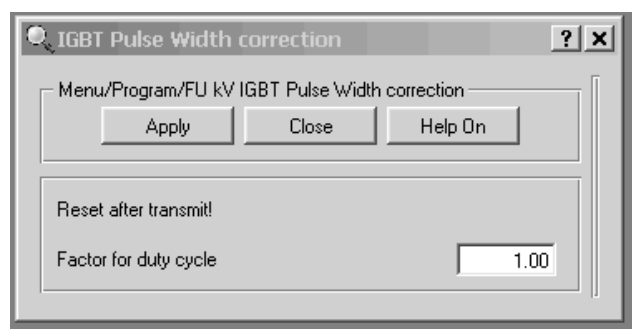



- If the  $U_{COMP}$  value does not match the requirements type in another factor duty cycle value.

Transmit the factor by clicking on “Apply” with the left mouse button and push the active RGDV button to get the new value validated.

- Switch an exposure.

The values are measured in the stationary condition. The transient behavior at the beginning of the exposure is not taken into account.

**Result:** In standby the  $U_{COMP}$  value is at about +11V, during exposure the mean value  $U_{COMP}$  must be as given in table 1, refer to figure 4:

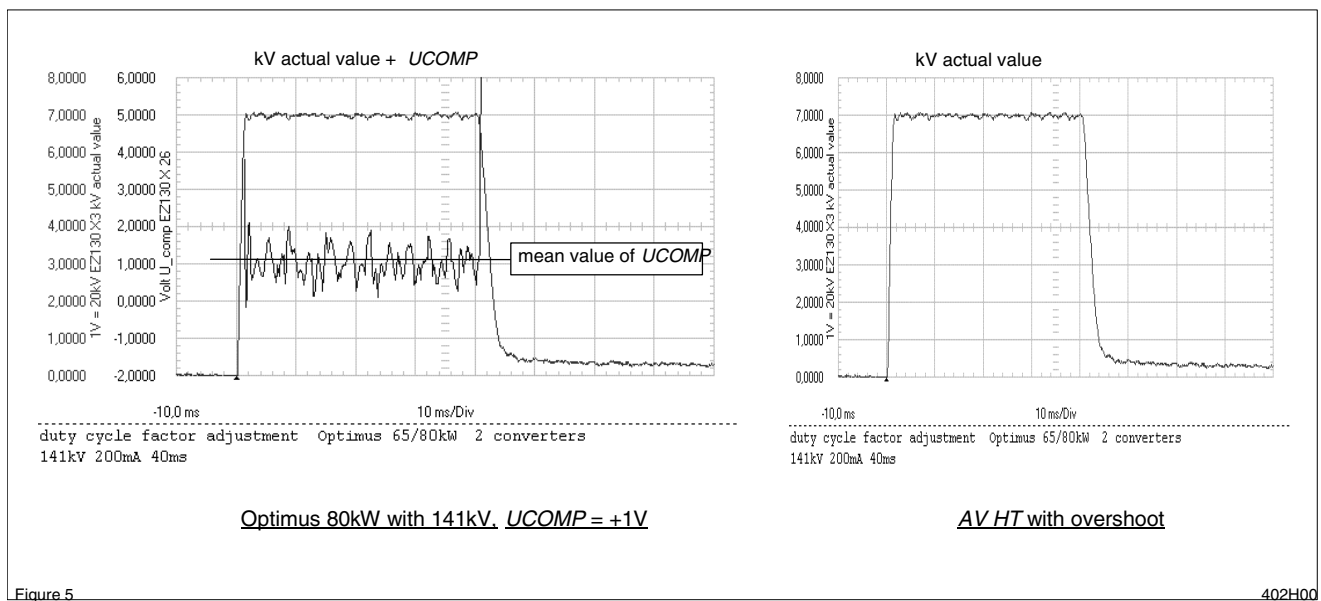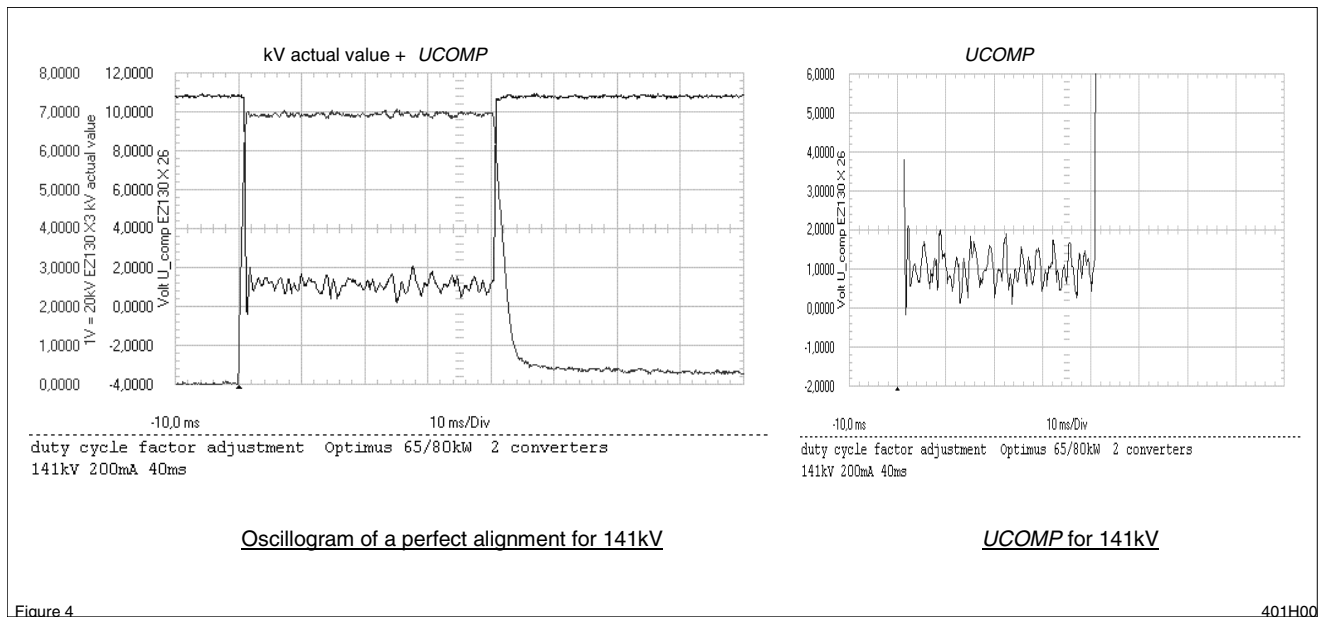

**125kV setting**

- Read the mean value of  $U_{COMP}$  for 125kV (in principle figure 4 or 5).
- Correct the factor duty cycle until  $U_{COMP}$  meets the required reference of 0V.

kV setpoint	mA setpoint	PCB type	$U_{COMP}$	Tolerance	kV peak of AV HT	Factor duty cycle:	Date
125kV	100mA	PCB kV_control 3:	+0V	+1V / -0.5V	125kV		

Table 1: Factor duty cycle, 125kV limit

Example how to correct the factor duty cycle:

- If the mean value of  $U_{COMP}$  is:
  - > +1V      **increase** the factor duty cycle in steps of 0.01
  - < -0.5V    **decrease** the factor duty cycle in steps of 0.01
- Check also the kV peak value AV HT (not the overshoot), it must be **125kV** for **125kV** setpoint.
- Turn OFF the generator and ENF1.
- Remove short link EZ130 X23 GAIN\_IN.
- Record the results in table 1

**2.4.3. Disposal**

- Dispose:
  1. PCB EZ 130 (4512 108 0908x)
  2. Old KV-Control RAD Software (4512 113 201xx)
  3. New KV-Control R/F Software (4512 113 26221) not used

according to local regulations.

**2.4.4. Final test**

- Switch ENF1 ON.
- Switch the generator ON.
- Test all applicable functions.

## 2.5. INSTALLATION IN OPTIMUS R/F OR OPTIMUS C GENERATORS

### 2.5.1. Installation of the R/F software

- Install the KV-Control R/F software 4512 113 26221 on the new PCB KV-EZ130 (4512 178 00261) in socket D800.
- Install the new PCB EZ130 (4512 178 00261) in the generator.

### 2.5.2. Adjustment of the factor for duty cycle

#### 2.5.2.1. General information


The actual value of the set kV must be attained at least after 2ms. During the kV rise phase there must be neither kV break-in nor a kV overshoot.

The factor duty cycle is based on an adapted tube and determines at local mains voltage and mains resistance conditions:

- kV rise phase
- kV behavior during the exposure in falling load technique.

The factor duty cycle is stored in the memory of PCB CU EZ139.


Refer to figure 1:

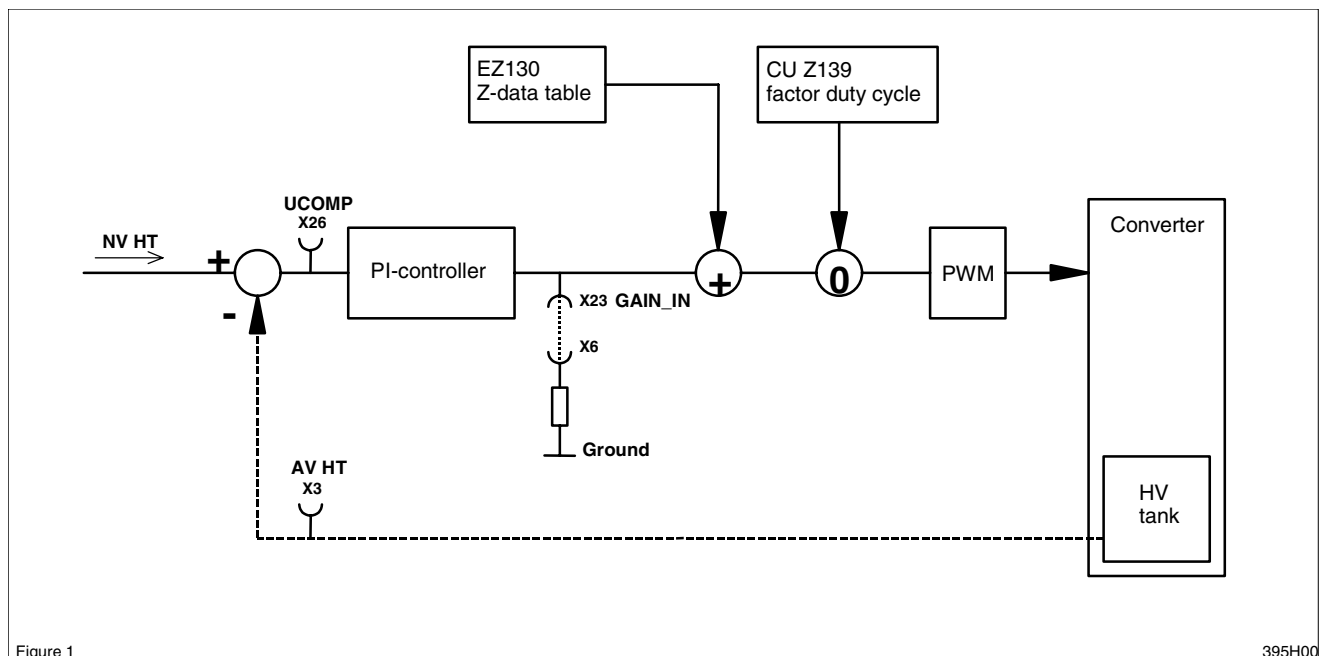

During alignment this factor duty cycle must be entered via AGenT. The influence of this factor as a correction value for the Z-data table is monitored as the  $U_{COMP}$  signal, since the PI-controller is deactivated by the grounded  $GAIN\_IN$  signal.

### 2.5.2.2. Connecting and setting the scope

For connections see figure 2:

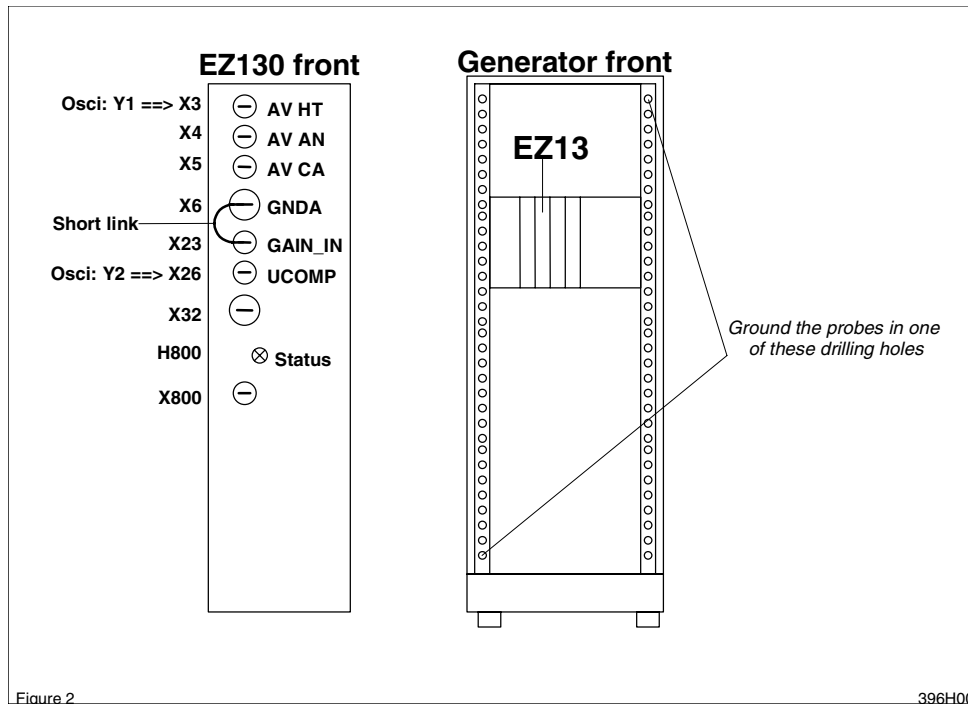

Channel 1 = EZ130 X3 ---> AV HT ---> 20kV/V ---> 1V/div --> Zero-line at bottom of screen  
 Probe GND = one of the drilling holes at the front cabinet chassis

Channel 2 = EZ130 X26 --->  $U_{COMP}$  ---> 1V/div ---> Zero-line 2 div from bottom of screen  
 Probe GND = one of the drilling holes at the front cabinet chassis

Trigger = external (preferred) ---> CTRL\_X\_C/ ---> backpanel EZX74 / negative slope  
 or = internal channel 1 ---> AV HT ---> EZ130 X3 / positive slope at +3V  
 Probe GND = one of the drilling holes at the front cabinet chassis

Time base = 5 or 10ms/div ---> trigger delay -1div



#### NOTE

*A digital scope should not have any other ground connection than the ground of the three probes at the drilling holes at the front generator chassis.*

*A mains-driven scope must be isolated from earth potential, otherwise it might display artefacts.*

### 2.5.2.3. Deactivating the kV controller

- Connect EZ130 X23 GAIN\_IN and X6 on board with a short link (use a short wire).



#### CAUTION

*This alignment requires exposures with high kV.  
Be sure the tube has been warmed up before.*

---

### 2.5.2.4. Setting of exposure data

- Switch ENF1 and the generator ON.
- Set 125kV.



#### NOTE

*Any tube arcing during this adjustment requires the execution of the tube conditioning next as described in chapter 3 (see page 22).*

*Disconnect the short link between X23 and X6.*


*Start over this adjustment from chapter 2.5.2.3. onwards if the tube conditioning was successful.*

---

- Set kV and mA values according to the programmed tube limits:  
**125kV:** 100mA at kV\_3 (50kW)
- Set the exposure time: 40ms



### 2.5.2.5. Adjustment of the factor for duty cycle

- Adjust the factor duty cycle via service software AGenT by measuring  $U_{COMP}$  with the scope.
- Connect the service PC and start Agent:  
Select menu:  
*Program / FU KV IGBT Pulse Width correction*
- Set the starting value factor duty cycle to **1.00**:

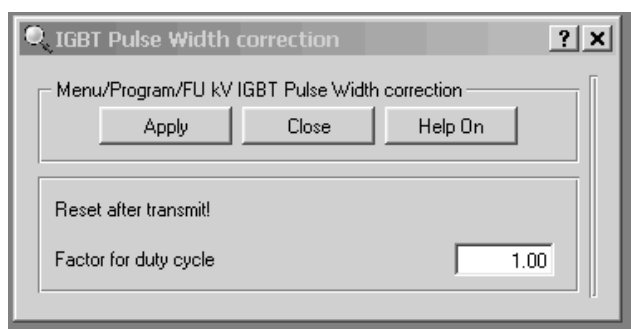

- If the  $U_{COMP}$  value does not match the requirements type in another factor duty cycle value.  
Transmit the factor by clicking on “Apply” with the left mouse button and push the active RGDV button to get the new value validated.
- Switch an exposure.  
The values are measured in the stationary condition. The transient behavior at the beginning of the exposure is not taken into account.

**Result:** In standby the  $U_{COMP}$  value is at about +11V, during exposure the mean value  $U_{COMP}$  must be as given in table 1, refer to figure 4:

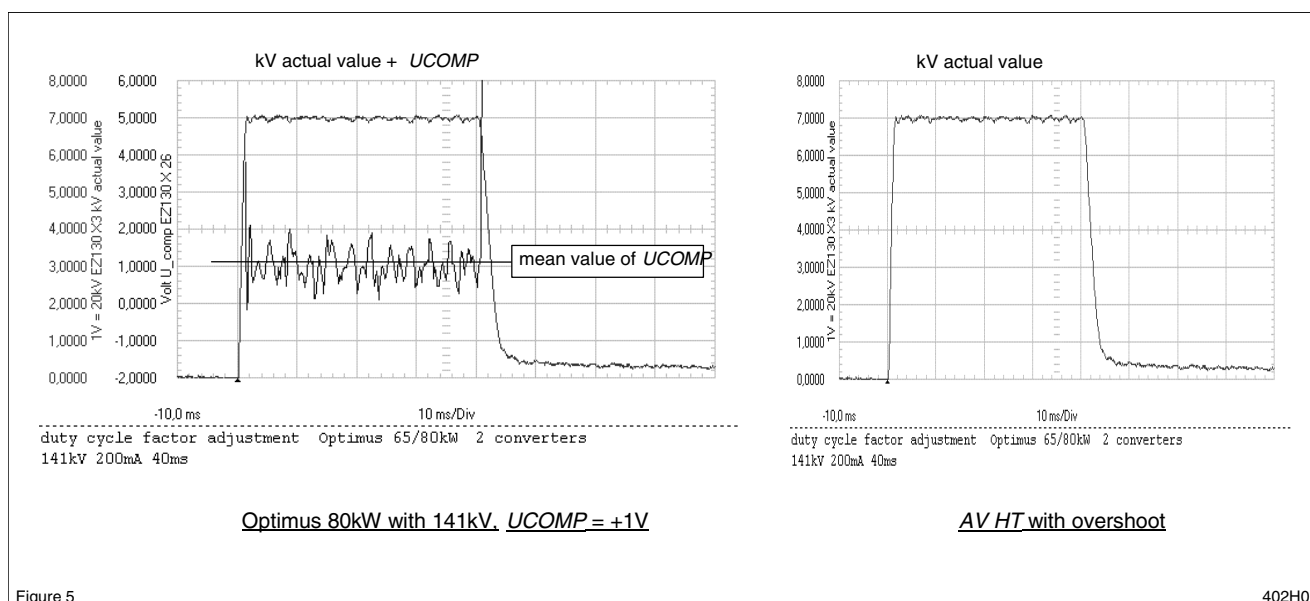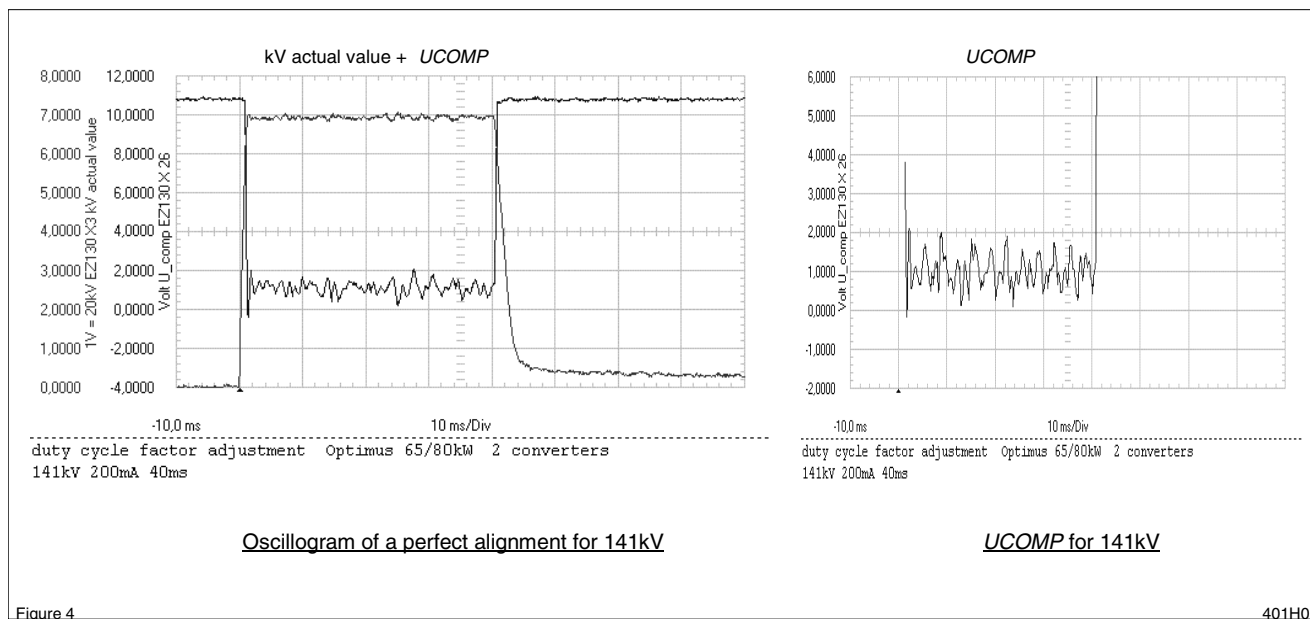



**125kV setting**

- Read the mean value of  $U_{COMP}$  for 125kV (in principle figure 4 or 5).
- Correct the factor duty cycle till  $U_{COMP}$  meets the required reference of 0V.

kV setpoint	mA setpoint	PCB type	$U_{COMP}$	Tolerance	kV peak of AV HT	Factor duty cycle:	Date
125kV	100mA	PCB kV_control 3:	+0V	+1V / -0,5V	125kV		

Table 1: Factor duty cycle, 125kV limit

- If the mean value of  $U_{COMP}$  is:
 

>	+1V	<b>increase</b> the factor duty cycle in steps of 0.01
<	-0.5V	<b>decrease</b> the factor duty cycle in steps of 0.01
- Check also the kV peak value *AV HT* (not the overshoot), it must be **125kV** for **125kV** setpoint.
- Turn OFF the generator and ENF1.
- Remove short link EZ130 X23 *GAIN\_IN*.
- Record the results in table 2.

**2.5.3. Disposal**

- Dispose:
  1. PCB EZ130 4512 108 0908x
  2. Old KV-Control software
  3. New KV-Control RAD software (4512 113 20141) not used

according to local regulations.

**2.5.4. Final test**

- Switch ENF1 ON.
- Switch the generator ON.
- Test all applicable functions.

## 2.6. INSTALLATION IN OPTIMUS RAD FOR US ARMY APPLICATION GENERATORS



### NOTE

*This OPTIMUS RAD generator is a special generator that is operated with an R/F high-voltage tank with expansion bellow (black). Therefore PCB EZ130 is **not** equipped with KV-Control RAD software, but with KV-Control R/F software.*

### 2.6.1. Installation of the R/F software

- Install the KV-Control R/F software 4512 113 26221 on the new PCB EZ130 (4512 178 00261) in socket D800.
- Install the new PCB EZ130 (4512 178 00261) in the generator.

### 2.6.2. Adjustment of the factor for duty cycle


#### 2.6.2.1. General information

The actual value of the set kV must be attained at least after 2ms. During the kV rise phase there must be neither kV break-in nor a kV overshoot.

The factor duty cycle is based on an adapted tube and determines at local mains voltage and mains resistance conditions:

- kV rise phase
- kV behavior during the exposure in falling load technique.


The factor duty cycle is stored in the memory of PCB CU EZ139. Refer to figure 1:

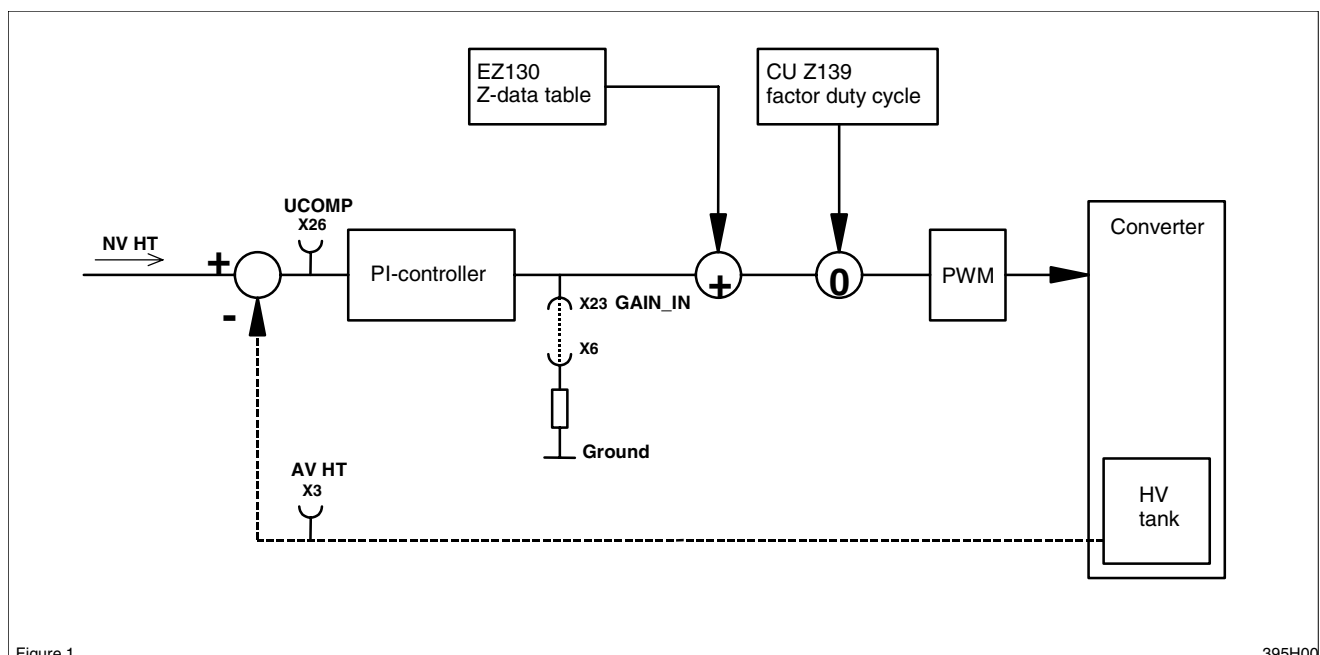

During alignment this factor duty cycle must be entered via AGenT. The influence of this factor as a correction value for the Z-data table is monitored as the  $U_{COMP}$  signal, since the PI-controller is deactivated by the grounded  $GAIN\_IN$  signal.

### 2.6.2.2. Connecting and setting the scope

For connections see figure 2:

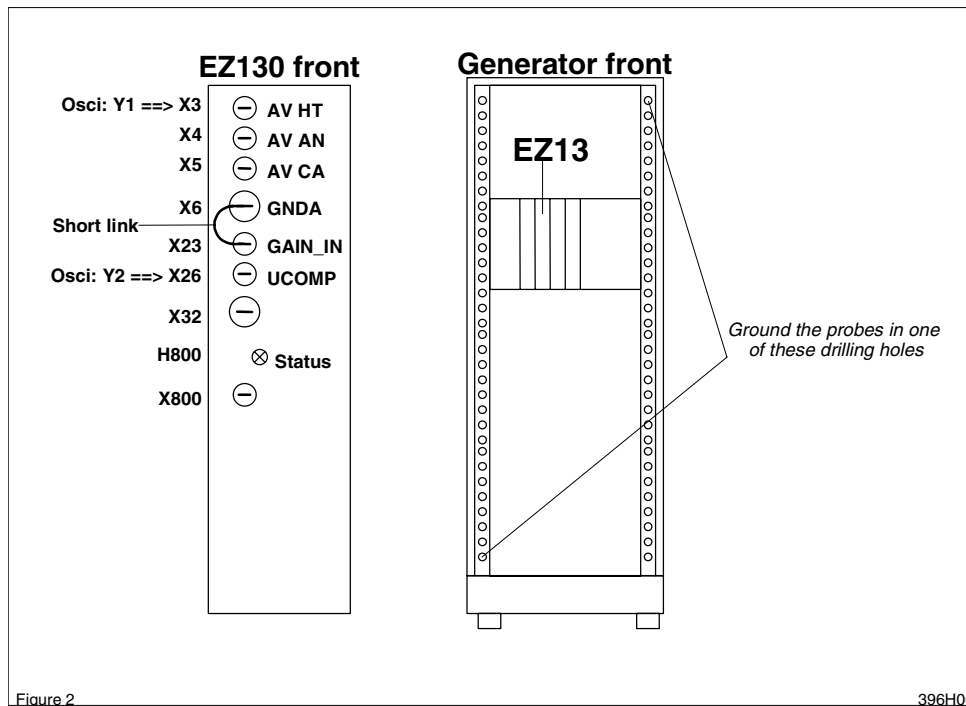

Channel 1 = EZ130 X3 ---> AV HT ---> 20kV/V ---> 1V/div --> Zero-line at bottom of screen

Probe GND = one of the drilling holes at the front cabinet chassis

Channel 2 = EZ130 X26 --->  $U_{COMP}$  ---> 1V/div ---> Zero-line 2 div from bottom of screen

Probe GND = one of the drilling holes at the front cabinet chassis

Trigger = external (preferred) ---> CTRL\_X\_C/ ---> backpanel EZX74 / negative slope

or = internal channel 1 ---> AV HT ---> EZ130 X3 / positive slope at +3V

Probe GND = one of the drilling holes at the front cabinet chassis

Time base = 5 or 10ms/div ---> trigger delay -1div



#### NOTE

*A digital scope should not have any other ground connection than the ground of the three probes at the drilling holes at the front generator chassis.*

*A mains-driven scope must be isolated from earth potential, otherwise it might display artefacts.*

### 2.6.2.3. Deactivating the kV controller

- Connect EZ130 X23 GAIN\_IN and X6 on board with a short link (use a short wire).



#### CAUTION

*This alignment requires exposures with high kV.  
Be sure the tube has been warmed up before.*

---

### 2.6.2.4. Setting of exposure data

- Switch ENF1 and the generator ON.
- Set 125KV.



#### NOTE

*Any tube arcing during this adjustment requires the execution of the tube conditioning next as described in chapter 3 (see page 22).*

*Disconnect the short link between X23 and X6.*


*Start over this adjustment from chapter 2.6.2.3. onwards if the tube conditioning was successful.*

---

- Set kV and mA values according to the programmed tube limits:  
**125kV:** 100mA at kV\_3 (50kW)
- Set the exposure time: 40ms

### 2.6.2.5. Adjustment of the factor for duty cycle

- Adjust the factor duty cycle via service software AGenT by measuring  $U_{COMP}$  with the scope.
- Connect the service PC and start AGenT:  
Select menu:  
*Program / FU KV IGBT Pulse Width correction*
- Set the starting value factor duty cycle to **1.00**:

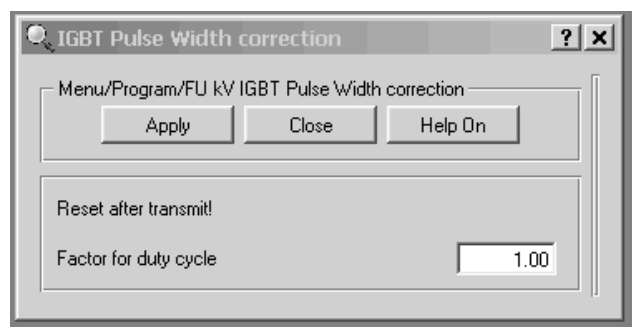

- If the  $U_{COMP}$  value does not match the requirements type in another factor duty cycle value.

Transmit the factor by clicking on “Apply” with the left mouse button and push the active RGDV button to get the new value validated.

- Switch an exposure.

The values are measured in the stationary condition. The transient behavior at the beginning of the exposure is not taken into account.

**Result:** In standby the  $U_{COMP}$  value is at about +11V, during exposure the mean value  $U_{COMP}$  must be as given in table 1, refer to figure 4:

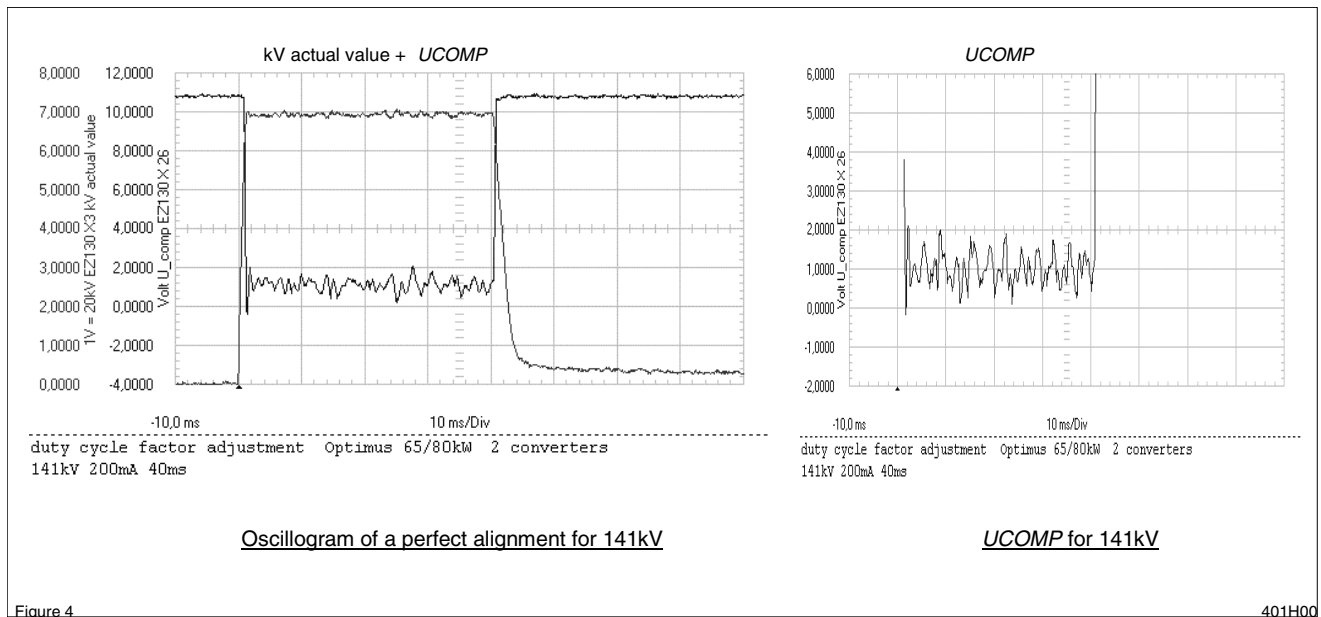


Figure 4

401H00

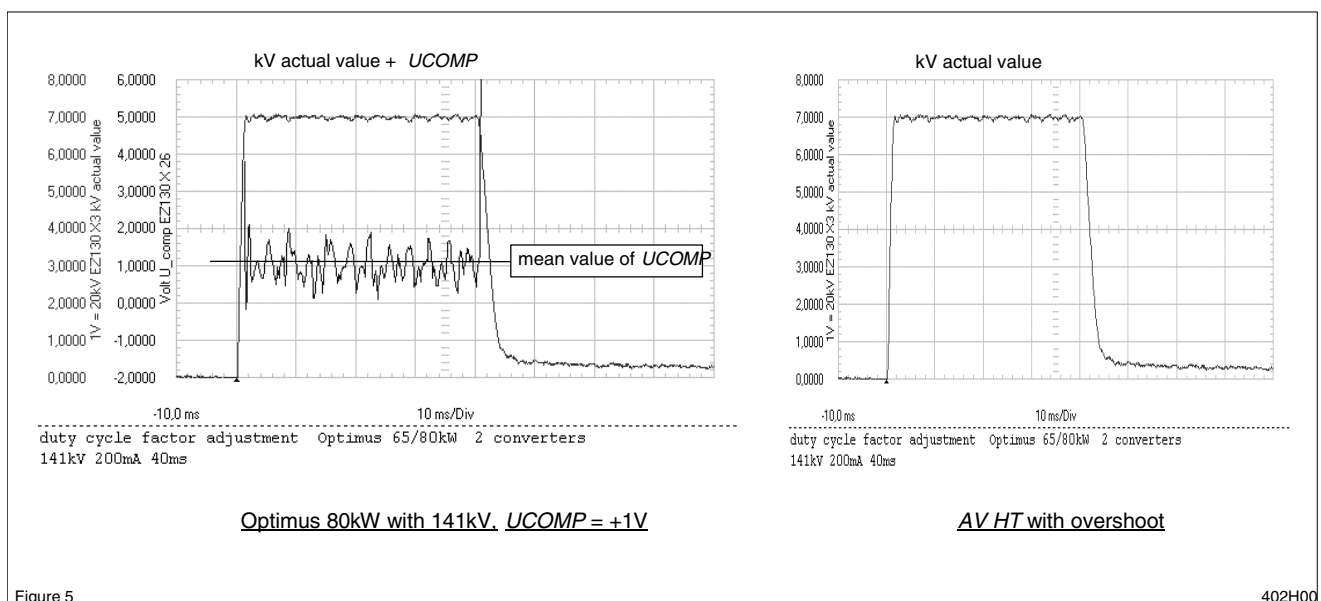


Figure 5

402H00

**125kV setting**

- Read the mean value of  $U_{COMP}$  for 125kV (in principle figure 4 or 5).
- Correct the factor duty cycle till  $U_{COMP}$  meets the required reference of 0V.

kV setpoint	mA setpoint	PCB type	$U_{COMP}$	Tolerance	kV peak of AV HT	Factor duty cycle:	Date
125kV	100mA	PCB kV_control 3:	+0V	+1V / -0,5V	125kV		

Table 1: Factor duty cycle, 125kV limit

Example how to correct the factor duty cycle:

- If the mean value of  $U_{COMP}$  is:
 

$> +1V$	<b>increase</b> the factor duty cycle in steps of 0.01
$< -0.5V$	<b>decrease</b> the factor duty cycle in steps of 0.01
- Check also the kV peak value *AV HT* (not the overshoot), it must be **125kV** for **125kV** setpoint.
- Turn OFF the generator and ENF1.
- Remove the short link EZ130 X23 *GAIN\_IN*.
- Record the results in table 1.

**2.6.3. Disposal**

- Dispose:
  1. EZ130 4512 108 09083
  2. Old KV-Control R/F software (4512 113 26212)
  3. New KV-Control RAD software (4512 113 20141) not used

according to local regulations.

**2.6.4. Modification of the documentation**

- Replace the pink page in the generator documentation (4512 984 28471) with the one supplied in this manual.

**2.6.5. Final test**

- Switch ENF1 ON.
- Switch the generator ON.
- Test all applicable functions.

### 3. TUBE CONDITIONING

#### 3.1. GENERAL



#### WARNING

*Radiation is released during the conditioning procedure!*

---

The generator must be in the **READY** state, i.e. the green ring at the desk must be illuminated!

#### 3.2. PRECONDITIONS / PROGRAM SETTINGS

- Switch OFF the generator.

Preparation of generators which are connected via a CAN interface:

- BuckyDiagnost TH and TH2
- DigitalDiagnost
- Thoravision
- EasyDiagnost with bucky unit

- Disconnect the following plugs.

System	Connector		
	EZX23 signal bus	EZX42 or EZX42-1 system CAN	EZX43 or EZX43-1 system CAN
BuckyDiagnost TH / TH2	X		X
DigitalDiagnost	X	X	X
Thoravision	X	X	X
EasyDiagnost with bucky unit	X	X	X

- Switch ON the generator.



#### NOTE

*The programming procedure must not be started before relay ENK1 has been energized at least 2 minutes after the generator has been switched on.*

---

- Perform the following program settings temporarily for each tube connected to one of the assigned RGDVs = Free cassette  
Select menu AGenT:  
*Program / RGDV set A + B / RGDV 1 ... 8 / Data Set A*

Program setting	Temporarily	Original tube
Enable handswitch .....	YES	
Syncmaster present	NO	
Exposure switch type	Double Step	
Exposure series / Tomo .....	YES	
Mounted radiographic .....	NONE	

- Reset the generator.
- Select appropriately programmed RGDV = "Free cassette" for the tube to be conditioned.

### 3.3. PROCEDURE

- Select **large** focus only.



#### NOTE

*The generator must be in the READY state.*

- Run reconditioning procedure for an adapted tube. Refer to the following table, left column "Tube adapted".
- or
- Run conditioning procedure for a new or non-adapted tube. Refer to the following table right column "Tube not adapted".
- It is recommended that the high voltage be monitored during conditioning.  
Connect the scope:  
Channel1: kV AV HT at EZ130 X3 (1V/div), scale: 20kV/V  
Trigger external: CTRL\_X\_C/ at backpanel EZ X74, negative slope  
Time base: 2ms/div

- In case of problems like tube arcing see the following flowchart EXPOSURE SEQUENCE as an example.  
The flowchart is valid for the applicable kV range only, e.g.:  
109kV is the max. kV value for normal application, set the next higher kV step = 117kV.



#### NOTE

*Refer to flowchart EXPOSURE SEQUENCE.*


*If the tube arcs at a certain kV value, switch another three exposures with same parameters and 10s pause between subsequent exposures. In case of success (no arcing anymore) continue with next kV step of the following table.*

*If the last exposure still arcs go one kV step back and follow the normal procedure. If this routine has been performed three times without improvement: ==> **Replace the tube!***



Exposure parameters for conditioning					
Tube adapted			# Exposures	Tube not adapted	
kV	mA	ms		kV	mAs
80	10	50	< 1 >	80	0.5
80	10	500	< 1 >	80	5
80	200	250	< 1 >	80	50
10 seconds pause				10 seconds pause	
80	max. mA	100	< 1 >	80	100
1 minute pause				1 minute pause	
90	10	50	< 1 >	90	0.5
90	10	500	< 1 >	90	5
90	200	250	< 1 >	90	50
10 seconds pause				10 seconds pause	
90	max. mA	100	< 1 >	90	100
1 minute pause				1 minute pause	
100	10	50	< 1 >	100	0.5
100	10	500	< 1 >	100	5
100	200	250	< 1 >	100	50
10 seconds pause				10 seconds pause	
100	max. mA	100	< 1 >	100	100
1 minute pause				1 minute pause	
110	10	50	< 1 >	110	0.5
110	10	500	< 1 >	110	5
110	200	250	< 1 >	110	50
10 seconds pause				10 seconds pause	
110	max. mA	100	< 1 >	110	100
1 minute pause				1 minute pause	
120	10	50	< 1 >	120	0.5
120	10	500	< 1 >	120	5
120	200	250	< 1 >	120	50
10 seconds pause				10 seconds pause	
120	max. mA	100	< 1 >	120	100
1 minute pause				1 minute pause	
130	10	50	< 1 >	130	0.5
130	10	500	< 1 >	130	5
130	200	250	< 1 >	130	50
10 seconds pause				10 seconds pause	
130	max. mA	100	< 1 >	130	100
1 minute pause				1 minute pause	

Exposure parameters for conditioning					
Tube adapted			# Exposures	Tube not adapted	
kV	mA	ms		kV	mAs
140	10	50	< 1 >	140	0.5
140	10	500	< 1 >	140	5
140	200	250	< 1 >	140	50
10 seconds pause				10 seconds pause	
140	max. mA	100	< 1 >	140	100
1 minute pause				1 minute pause	
145	10	50	< 1 >	145	0.5
145	10	500	< 1 >	145	5
145	200	250	< 1 >	145	50
10 seconds pause				10 seconds pause	
145	max. mA	100	< 1 >	145	100
1 minute pause				1 minute pause	
148	10	50	< 1 >	148	0.5
148	10	500	< 1 >	148	5
148	200	250	< 1 >	148	50
10 seconds pause				10 seconds pause	
148	max. mA	100	< 1 >	148	100
1 minute pause				1 minute pause	
150	10	50	< 1 >	150	0.5
150	10	500	< 1 >	150	5
150	200	250	< 1 >	150	50
10 seconds pause				10 seconds pause	
150	max. mA	100	< 1 >	150	100
1 minute pause				1 minute pause	

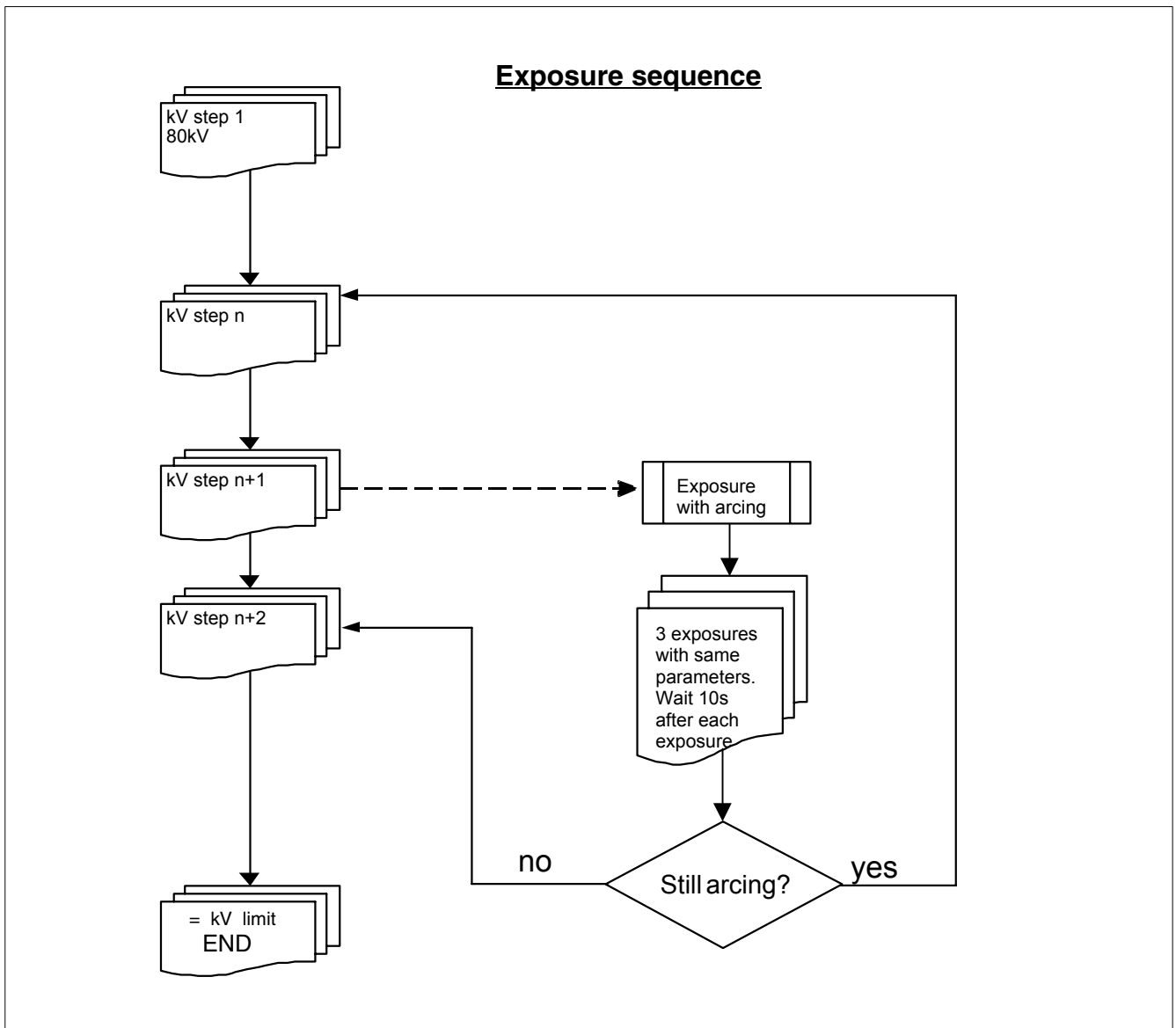**NOTE**

*If a tube arcs at any kV value which is not required for application, program the max. application kV value with AGenT:*

*Program / Tubes / Tube Limits / Max. Tube Voltage Limit [kV] / [117]*

As the max. kV value has decreased now, the field ADAPTED TO [kV] displays the max. value after adaptation as well.

- Set the RGDV programming to the original status if no adaptation procedure has to be executed.
- Reset the generator.